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Intended use


Blood Gas Monovette® systems are blood collection systems and are used together with needles, adapters or catheters for arterial and in rare cases venous blood collection for blood gas analysis. Blood Gas Monovette® systems are used to collect blood samples and for point-of-collection testing (POCT), as well as for direct determination of blood gases, pH or electrolytes on a blood gas analyser within 15 minutes* after the blood has been collected. The products are intended for use in a professional environment by qualified medical and laboratory personnel.

Product description

Blood Gas Monovette® systems consist of a plastic container, a plunger, a male Luer connector (screw cap), an orange plastic cap and calcium-balanced heparin (50 IU/ml blood, spray-dosed in droplet form**) as a preparation. The blood collection system is available in a sterile inside version or in an individually packaged version sterile on the inside and outside.

The accessory is a venting device (art. no. 14.1148), which like the membrane adapter (art. no. 14.1112) is available separately, but also pre-assembled on the Blood Gas Monovette® system.

NOTE: No lithium determinations may be carried out from the Blood Gas Monovette®.
* Thomas Clinical Laboratory Diagnostics: www.clinical-laboratory-diagnostics-2020.com/chapter-9-Acid-base-balance-and-blood-gases/chapter-9.5-Specimen
** Gruber et al. Heparin release is insufficient in syringes with platelets as heparin source CCA 2008; doi:10.1016/j.cca.2008.04.029.

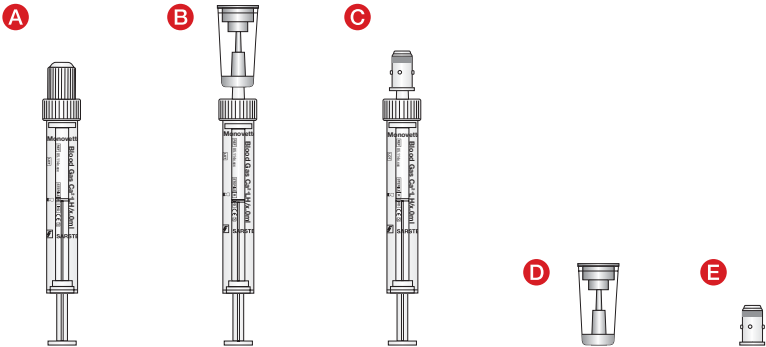Type	Description
A	Blood Gas Monovette® with cap, sterile
B	Blood Gas Monovette® with venting device, sterile
C	Blood Gas Monovette® with membrane adapter, sterile
D	Venting device, non-sterile
E	Membrane adapter, sterile

Safety information and warnings

- General precautions: Use gloves and other general personal protective equipment to protect yourself against blood and possible exposure to pathogens transmitted by blood.
- Handle all biological specimens and sharp blood collection utensils (needles) according to the guidelines and procedures of your facility. In case of direct contact with biological specimens or a needlestick injury, consult a doctor due to the risk of transmission of HIV, HCV, or HBV or other infectious diseases. Use safety needles that have a built-in needle guard. Always follow the safety guidelines and procedures of your facility.
- If blood is withdrawn via an arterial or intravenous (IV) catheter, it must be ensured that the catheter is properly flushed (= cleaned from the IV solution) according to the rules of the facility before the blood collection is started. Properly flushing the line prevents faulty laboratory results.
- Underfilling the Blood Gas Monovette® leads to an incorrect ratio of blood to preparation/additive and can lead to incorrect analysis results. Overfilling is not possible due to the shape.
- The blood removed and processed with the Blood Gas Monovette® is not intended for re-injection into the human body.
- Dispose of all sharp objects (e.g. needles) for blood collection in suitable sharps containers.
- Do not use the product after the expiry date. The expiry date of the product is the last day of the month and year indicated.

Storage

The Blood Gas Monovette® should be stored at room temperature.

Transport

The product serves as a primary container for in-house, timely transport for point-of-collection blood gas analysis at room temperature.

Limitations

Storing the sample for blood gas analysis is not possible, since the analysis takes place near the point of collection and should be finished after 15 minutes. Intended only for in-house transport in a timely manner.

NOTE: No lithium determinations may be carried out from the Blood Gas Monovette®.

Collecting and handling specimens

READ THIS DOCUMENT COMPLETELY BEFORE STARTING A VENIPUNCTURE.

Equipment needed to collect a specimen:

1. Gloves, gown, goggles or other suitable protective wear to protect against bloodborne pathogens or potentially infectious materials.
2. Prepare the Blood Gas Monovette® and venting device.
3. Labels for specimen identification.
4. S-Monovette® safety needle and membrane adapter, arterial catheter or intravenous catheter.
5. Disinfectant for cleaning the puncture site (follow the facility's guidelines for preparing the puncture site for taking blood samples).
6. Dry, clean swabs.
7. When collecting the sample from the catheter: Stoppers or similar for catheters according to your facility's guidelines.
8. Vein tourniquet for venous blood collection.
9. Plaster or dressing.
10. Sharps disposal container for the safe disposal of used material.

Preventing reflux

The Blood Gas Monovette® contains a chemical additive. To prevent reflux, the following should be noted:

1. Use only the materials recommended in the instructions for use.
2. Position the patient so that there is easy access to the artery or to the vein and the arm or other puncture site on the patient is directed downwards as much as possible.
3. In addition, keep the Blood Gas Monovette® below the arm level/ puncture range as far as possible.
4. It must be ensured that the contents of the Blood Gas Monovette® do not touch the Luer connection or the screw cap during puncture.
5. For venous blood collection, loosen the vein tourniquet while the blood is flowing into the Blood Gas Monovette®.

Vascular puncture / vascular access


General instructions:

WEAR GLOVES TO MINIMISE THE RISK OF EXPOSURE DURING BLOOD COLLECTION.

A: Arterial puncture / connection with arterial catheter

Select arterial puncture site / arterial catheter and prepare according to the facility's guidelines.


1. Pull off the orange cap and set aside.
- 2a. Grip the Blood Gas Monovette® on the container body (NOT on the screw cap) and connect to the membrane adapter by pushing the Blood Gas Monovette® into the female Luer connection of the membrane adapter and securing it with a clockwise movement. This system is completed with a SARSTEDT needle (e.g. the S-Monovette® (safety) needle), or
- 2b. connect the Blood Gas Monovette® to the female Luer connection of an arterial catheter by pushing the Blood Gas Monovette® into the female Luer connection of the catheter and securing it in place with a clockwise movement.
- 3a. Puncture the artery at an angle of 45°. Subsequently, draw the blood by gently pulling the plunger of the Blood Gas Monovette® until it stops. When the blood flow stops, the Blood Gas Monovette® is disconnected and the puncture site treated in accordance with the facility's guidelines, or
- 3b. the blood is drawn by gently pulling the plunger of the Blood Gas Monovette® up to the stop. When the blood flow stops, disconnect the Blood Gas Monovette® and handle the catheter according to the facility's guidelines.
4. Immediately after blood collection, the air above the blood sample is removed with the attached venting device in order to avoid incorrect measurements due to air contamination.
5. After venting the sample, mix the blood sample directly with the closure cap** placed on it, by rolling the Blood Gas Monovette® back and forth between the palms of the hands 10x.
DO NOT mix by swinging overhead!
** The variant of the Blood Gas Monovette® with venting device (without cap) is mixed with the venting device attached.
6. Blood gas analyses should be performed up to 15 minutes after blood collection.

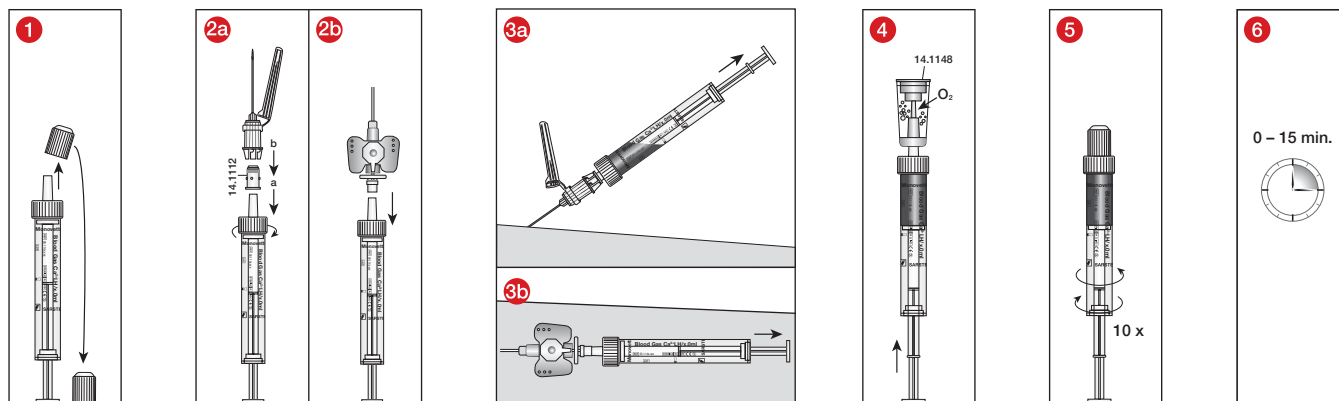B: Venipuncture /connection with intravenous catheter

Select venipuncture site / intravenous catheter and prepare according to the facility's guidelines.

- 1** See under A. Arterial puncture/Connecting with an arterial catheter.
- 2a** See A. Arterial puncture/connecting with the arterial catheter, or
- 2b** See A. Arterial puncture/connecting with the arterial catheter, but also with an intravenous catheter.
- 3a** Puncture the vein at a flat angle and then draw the blood by gently pulling back the plunger of the Blood Gas Monovette® until the stop.
When the blood flow stops, the Blood Gas Monovette® is disconnected and the puncture site treated in accordance with the facility's guidelines, or
- 3b - 6** see under A. Arterial puncture/Connecting with an arterial catheter.

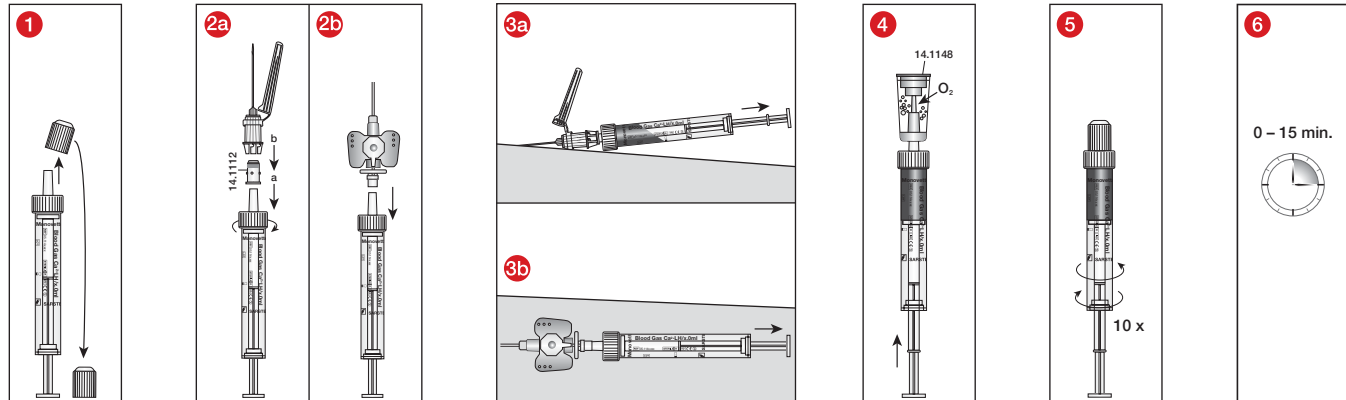Disposal

1. The general hygiene guidelines and regulations for the proper disposal of infectious material must be observed and complied with.
2. Disposable gloves prevent the risk of infection.
3. Contaminated or filled blood collection systems must be disposed of in suitable containers for hazardous biological waste that can be subsequently autoclaved and burned.
4. They must be disposed of in a suitable incinerator or by autoclaving (steam sterilisation).

Standards and recommendations

CLSI* C46 Blood Gas and pH Analysis and Related Measurements; Approved Standard

WHO/DIL/LAB/99.1 Rev02 "WORLD HEALTH ORGANIZATION, et al. Use of anticoagulants in diagnostic laboratory investigations. Geneva: World Health Organization, 2002".

*CLSI: Clinical and Laboratory Standards Institute

Key for symbols and labels:

	Article number
	Batch number
	Use by
	CE marking
	<i>In vitro</i> diagnostic device
	Follow the instructions for use
	If reused: Risk of contamination
	Keep away from sunlight
	Store in a dry place
	Manufacturer
	Country of manufacture

The following also applies to sterile products:

	Do not use if package is damaged
	Single sterile barrier system with external protective packaging
	Sterilisation by irradiation
or:	
	Sterile fluid path
	Do not resterilise

Technical modifications reserved.

All serious incidents relating to the product must be reported to the manufacturer and the competent national authority.

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