

# Urine Drainage

Product solutions for the individual patient care



# A06 W

- ✓ Needle-free and flat sampling site
- ✓ Attachment of the exchange bag with Bayonet-lock
- ✓ Universal attachment with elastic straps and additional hooks
- ✓ Exchange bags are available with 2,000 ml and 3,000 ml capacity




## Precise Reading



The precise and easily readable graduations facilitate accurate output readings. The total measuring volume of 500 ml and the clear and concise graduations enable an accurate balancing of large urine volumes. Volumes exceeding 500 ml are automatically diverted through the integral overflow chamber and drained into the collection bag.

## Hygienic & user-friendly exchange bag

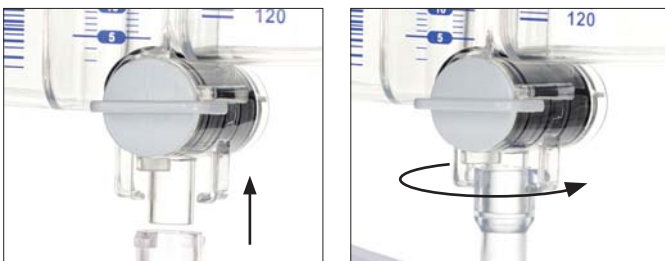

The Bayonet-Lock enables a quick and hygienic exchange of the bag and a secure connection to the measuring chamber.

DEHP free



## Needle-free sampling



Needle-free, flat sampling site for secure disinfection and sampling.

# UD 3PLUS

- ✓ Needle-free and flat sampling site
- ✓ Large double hook
- ✓ 2 tube lengths: 100 cm and 120 cm
- ✓ 2,000 ml capacity



## The System for the immobile patient

Closed drainage system for the hygienic care of catheterised patients with long-term confinement.

## Large double hook



The large double hook allows easy and secure attachment to the bed rails.

## Drip Chamber



With hydrophobic ventilation for optimal pressure compensation and urine flow. The drip chamber enables optimal protection against ascending migration of bacteria. The horizontal membrane reflux valve ensures optimal urine flow from small to large volumes and prevents urine reflux regardless of the drainage bag position.

DEHP free



## Needle-free sampling



Needle-free, flat sampling site for secure disinfection and sampling.

# TUR-BAG

- ✓ Needle-free and flat sampling site
- ✓ Available with and without pumpball
- ✓ 4,000 ml capacity



## The bag for bladder irrigation

A closed urine drainage system for the post-operative bladder irrigation.

## Pumpball



With pump ball for controlled postoperative bladder irrigation, especially after prostate resection.

## Needle-free sampling



Needle-free, flat sampling site for secure disinfection and sampling.

DEHP free



# DRAINAGE BAG

- ✓ Universal step connector with integrated pierceable membrane for urine sampling
- ✓ Large diameter drainage tubing for optimal urine flow
- ✓ 2,000 ml capacity



DEHP free

## Optimal Patientcare with the drainage bag

The drainage bag features an extra strong foil and a bacteria-proof, hydrophobic bag ventilation.

## T-Tap drainage valve



Easy-to-handle T-Tap drainage valve and effective protection from contamination.

## Needle-free sampling



Needle-free, flat sampling site for secure disinfection and sampling.




# URINE- AND SECRETIONBAGS

- ✓ universal step connector
- ✓ The urine- and secretion bags are available with 1,500 ml and 2,000 ml capacity
- ✓ Tube length 90 cm
- ✓ Available in sterile and non-sterile



## Anti-Reflux Valve

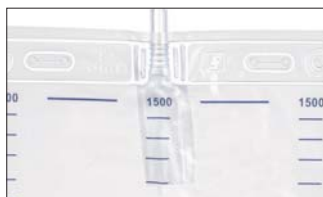

The anti-reflux valve reliably prevents urine from flowing back.

## Push-pull drainage valve




The push-pull drain valve enables quick and easy / hygienic emptying of the bag.

## Luer-Lock connector



The (male) Luer-Lock connector allows safe connection to (Gall)-drainages.

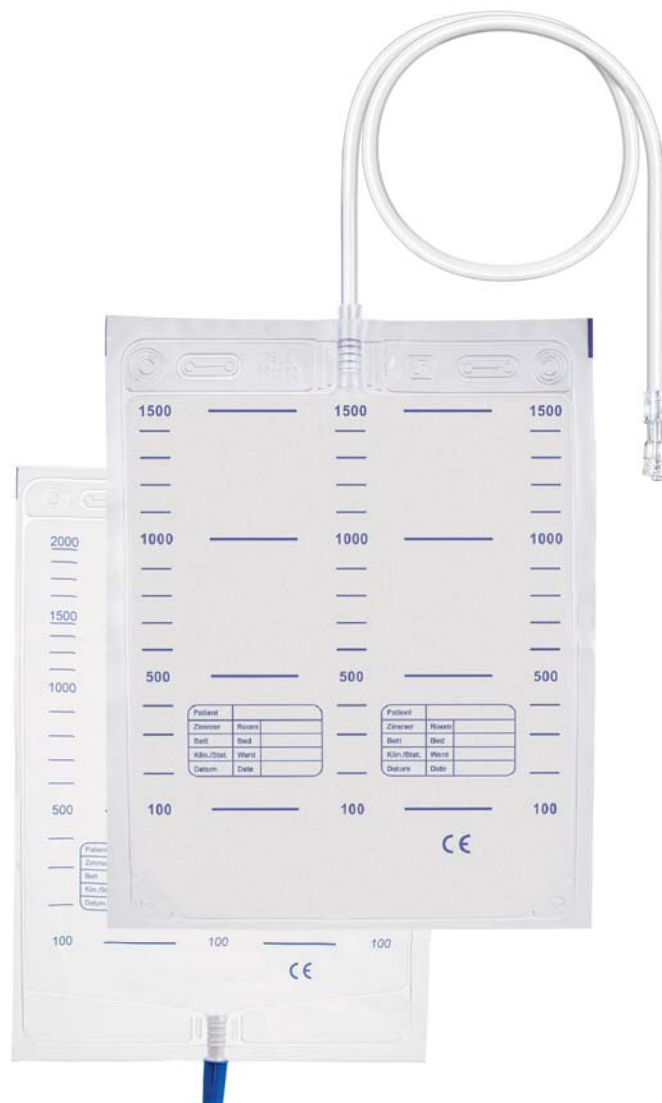

# UROSTRETCH

- ✓ No contact of the leg bag with the skin
- ✓ Full length fixation of leg bag surface
- ✓ Reduction of liquid noise



## UroStretch leg bag attachment



The stretch fabric, tubular bandage for the attachment of urine leg-bags. Usable with all leg bags up to 750 ml. UroStretch is reusable, washable in warm water (40°C), and offers high user comfort.



# ORDER INFORMATION



## Urine Measuring System A06 W

| Ref         | Description  | Pcs./Box |
|-------------|--|----------|
| 74.5681.200 | Urine measuring system A06 W 500 ml with 2.0 l exchange bag without drainage tap, with elastic straps and additional hooks, individually packed, sterile | 10       |
| 74.5684.001 | Exchange bag 2.0 l for A06 W without drainage tap, individually packed, sterile  | 50       |
| 74.5684.002 | Exchange bag 2.0 l for A06 W with drainage tap, individually packed, sterile   | 50       |
| 74.5685.001 | Exchange bag 3.0 l for A06 W without drainage tap, individually packed, sterile  | 50       |
| 74.5685.002 | Exchange bag 3.0 l for A06 W with drainage tap, individually packed, sterile   | 50       |

## UD 3plus

| Ref         | Description  | Pcs./Box |
|-------------|--|----------|
| 74.5330.500 | UD 3plus urine drainage system 2.0 l, with double hook, tube length 100 cm, individually packed, sterile | 20       |
| 74.5330.520 | UD 3plus urine drainage system 2.0 l, with double hook, tube length 120 cm, individually packed, sterile | 20       |

## TUR-BAG

| Ref         | Description   | Pcs./Box |
|-------------|---|----------|
| 74.5340.500 | TUR-BAG 4.0 l with pump ball, tube length 100 cm, individually packed, sterile    | 10       |
| 74.5340.501 | TUR-BAG 4.0 l without pump ball, tube length 100 cm, individually packed, sterile | 20       |

## Drainage bag

| Ref         | Description  | Pcs./Box |
|-------------|--|----------|
| 74.5485.005 | Drainage bag 2.0 l, with anti-reflux valve, T-Tap drainage valve, tube length 100 cm, individually packed, sterile | 50       |





### Urine- and secretionbags

| Ref         | Description  | Pcs. / Box |
|-------------|--|------------|
| 74.5218.101 | Urine bag 1.5 l, with anti-reflux valve, push-pull drainage valve, tube length 90 cm, non-sterile                  | 150        |
| 74.5228.101 | Urine bag 2.0 l, with anti-reflux valve, push-pull drainage valve, tube length 90 cm, non-sterile                  | 150        |
| 74.5228.005 | Urine bag 2.0 l, with anti-reflux valve, push-pull drainage valve, tube length 90 cm, individually packed, sterile | 100        |
| 74.5283.101 | Urine bag 1.5 l, with anti-reflux valve, tube length 90 cm, non-sterile  | 200        |
| 74.5283.005 | Urine bag 1.5 l, with anti-reflux valve, tube length 90 cm, individually packed, sterile                           | 100        |
| 74.5284.101 | Urine bag 2.0 l, with anti-reflux valve, tube length 90 cm, non-sterile  | 200        |
| 74.5284.005 | Urine bag 2.0 l, with anti-reflux valve, tube length 90 cm, individually packed, sterile                           | 100        |
| 74.5281.101 | Urine- and Secretion bag 1.5 l, tube length 90 cm, non sterile   | 200        |
| 74.5281.005 | Urine- and Secretion bag 1.5 l, tube length 90 cm, individually packed, sterile                                    | 100        |
| 74.5282.101 | Urine- and Secretion bag 2 l, tube length 90 cm, non-sterile   | 200        |
| 74.5282.005 | Urine- and Secretion bag 2 l, tube length 90 cm, individually packed, sterile                                      | 100        |
| 74.5286.005 | Gall secretion bag 1.5 l without drainage tap, with Luer-Lock, tube length 90 cm, individually packed, sterile     | 100        |

### UroStretch

| Ref         | Description                   | Pcs. / Box |
|-------------|-------------------------------|------------|
| 74.5767.001 | UroStretch, size: small       | 10         |
| 74.5767.002 | UroStretch, size: medium      | 10         |
| 74.5767.003 | UroStretch, size: large       | 10         |
| 74.5767.004 | UroStretch, size: extra large | 10         |

## SARSTEDT AG & Co. KG

P.O. Box 12 20

D-51582 Nümbrecht

Tel.: +49 2293 305 0

Fax: +49 2293 305 3992

[export@sarstedt.com](mailto:export@sarstedt.com)

[www.sarstedt.com](http://www.sarstedt.com)



**SARSTEDT**