

Sarpette[®] M

Systematic pipetting



Try it out for
free & without
obligation!



shop.sarstedt.us/sarpette



SARSTEDT

A PERFECT MATCH

- ✓ Compatible and compliant
- ✓ Robust and reliable
- ✓ Single source




Sarpette® M - Systematic Pipetting

SARSTEDT has been producing precision pipette tips for more than 50 years. As a result of our years of experience and technical expertise, and in the interest of being a single source supplier, we have developed the new Sarpette® M, an innovative microliter pipette that is perfectly tailored to today's applications in the laboratory.

Sustainability and environmental responsibility are increasingly important.


Adjustable tip ejector

Precision fit cone

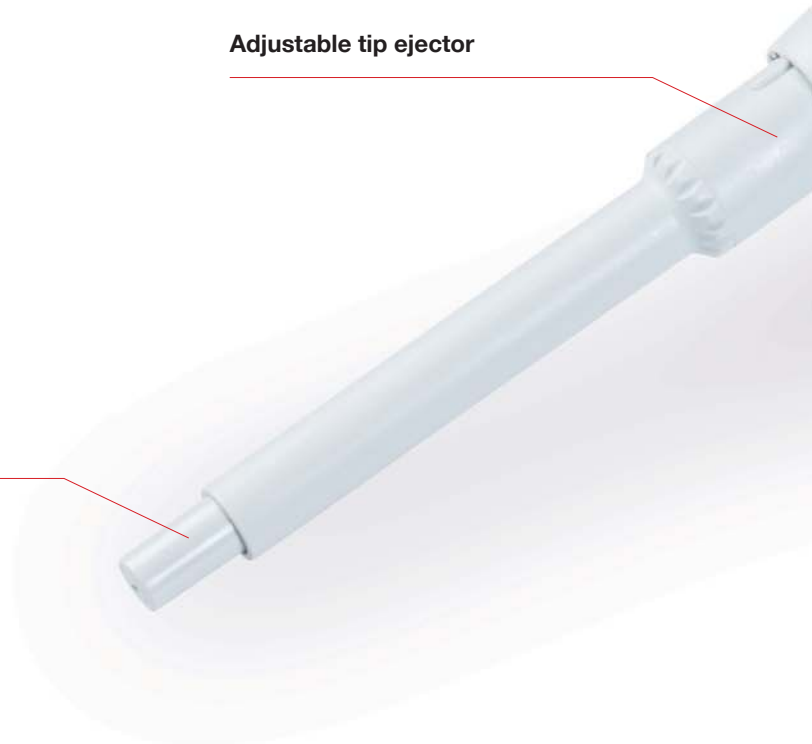

Sarpette[®] M

For high-precision pipetting

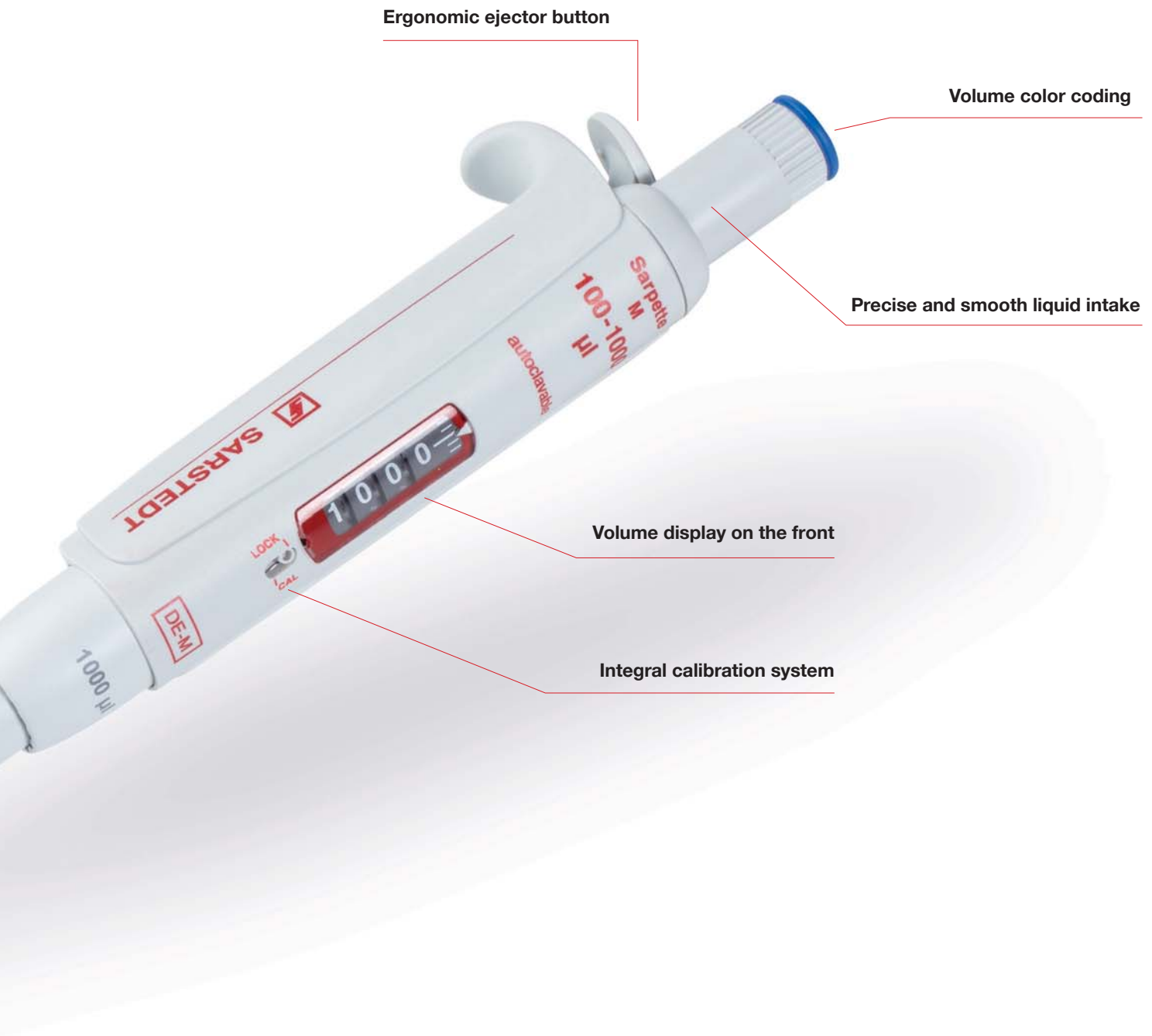

ERGONOMIC DESIGN

- ✓ Smooth and ergonomic pipetting button
- ✓ Convenient and comfortable design
- ✓ Flexible and adjustable tip ejection



User friendly design

Particularly ergonomic design

For the daily requirements in research, diagnostics, and any laboratory, SARSTEDT precision pipettes stand out for their excellent characteristics. Benefit from their exceptional lightweight design and durability.





Unique ejector position

The ejector directly above the hand rest enables a particularly comfortable user experience. Thanks to the large ejector push button and the smooth-running pipetting piston, the Sarpette® M guarantees ergonomic pipetting. Its light weight is an added benefit to prevent RSI (repetitive strain injury).



Ideal pipetting

Thanks to the flexibly adjustable tip ejector, you can easily adjust the ejector to the individual height of your pipette tips with just one turn and little effort. This makes the Sarpette® M compatible with a wide range of long and short pipette tips from multiple manufacturers.



RELIABLE

- ✓ Fast and accurate volume adjustment
- ✓ Accurate and precise pipetting
- ✓ Simple and safe calibration



Fast and safe

Easy volume adjustment

To set the desired volume, simply turn the free-moving and precise adjustment knob.

Audible click-stops and the rotatable thumb rest prevent unintended adjustment.





Precise and accurate

Precise pipetting results are a given with the Sarpette® M. SARSTEDT pipette tips and your Sarpette® M are a perfect match for all areas of application.



Calibration in a flash

Awkward calibration is a thing of the past. With the SARSTEDT Sarpette® M, you can now calibrate your pipettor yourself - quickly and safely. The precise system conforms to all international standards.



ROBUST

- ✓ Strong and durable construction
- ✓ Autoclavable and UV-resistant
- ✓ Convenient and quick cleaning



Only the best materials

Quick and easy disassembly

Disassemble your Sarpette® M within seconds. This allows you to clean and disinfect the pipette at any time. Convenient and quick cleaning makes contamination-free and safe pipetting even easier.





Fully autoclavable

The Sarpette® M can be easily autoclaved at 121°C. The pipette is extremely resistant to chemicals and UV exposure.



Stable and durable construction

Thanks to the premium materials and excellent workmanship of the Sarpette® M, we can offer a warranty period of 3 years.



ORDERING INFORMATION



Sarpette® M performance values

Order no.	Volume range	Graduation	Inaccuracy (E%)			Imprecision (CV%)		
			Min. vol.	Med. vol.	Max. vol.	Min. vol.	Med. vol.	Max. vol.
90.3100.002	0.1 - 2 µl	0.002 µl	<+/- 6.0 % ¹	<+/- 4.0 %	<+/- 2.0 %	< 5.0 % ¹	< 3.3 %	< 1.5 %
90.3100.010	0.5 - 10 µl	0.01 µl	<+/- 2.5 % ²	<+/- 1.8 %	<+/- 1.0 %	< 1.8 % ²	< 1.2 %	< 0.5 %
90.3100.020	2 - 20 µl	0.02 µl	<+/- 2.5 %	<+/- 1.8 %	<+/- 1.0 %	< 1.7 %	< 1.0 %	< 0.5 %
90.3100.100	10 - 100 µl	0.1 µl	<+/- 1.5 %	<+/- 1.2 %	<+/- 0.8 %	< 1.0 %	< 0.6 %	< 0.2 %
90.3100.200	20 - 200 µl	0.2 µl	<+/- 1.5 %	<+/- 1.1 %	<+/- 0.6 %	< 0.6 %	< 0.4 %	< 0.2 %
90.3100.000	100 - 1,000 µl	1 µl	<+/- 1.5 %	<+/- 1.0 %	<+/- 0.5 %	< 0.5 %	< 0.4 %	< 0.2 %
90.3100.555	0.5 - 5 ml	0.01 µl	<+/- 1.5 %	<+/- 1.1 %	<+/- 0.6 %	< 0.6 %	< 0.5 %	< 0.3 %
90.3100.111	1 - 10 ml	0.01 µl	<+/- 1.5 %	<+/- 0.7 %	<+/- 0.5 %	< 0.5 %	< 0.3 %	< 0.2 %

The performance values were determined using distilled water at a constant temperature (± 0.5 °C) between 20 and 25 °C in accordance with ISO 8655.

1) measured at 0.5 µl

2) measured at 1 µl

Sarpette® M8 / M12 performance values

Order no.	Volume		Number of channels	Inaccuracy (E%)			Imprecision (CV%)		
	Volume range	Graduation		Min. vol.	Med. vol.	Max. vol.	Min. vol.	Med. vol.	Max. vol.
90.3108.010	0.5 - 10 µl	0.01 µl	8	<+/- 3.5 % ¹	<+/- 2.5 %	<+/- 1.5 %	< 3.0 % ¹	< 2.0 %	< 1.0 %
90.3108.200	20 - 200 µl	0.2 µl	8	<+/- 0.9 % ¹	<+/- 0.8 %	<+/- 0.7 %	< 0.6 % ¹	< 0.5 %	< 0.3 %
90.3108.300	40 - 300 µl	0.4 µl	8	<+/- 1.0 % ¹	<+/- 0.9 %	<+/- 0.8 %	< 0.6 % ¹	< 0.5 %	< 0.3 %
90.3112.010	0.5 - 10 µl	0.01 µl	12	<+/- 3.5 % ¹	<+/- 2.5 %	<+/- 1.5 %	< 3.0 % ¹	< 2.0 %	< 1.0 %
90.3112.200	20 - 200 µl	0.2 µl	12	<+/- 0.9 % ¹	<+/- 0.8 %	<+/- 0.7 %	< 0.6 % ¹	< 0.5 %	< 0.3 %
90.3112.300	40 - 300 µl	0.4 µl	12	<+/- 1.0 % ¹	<+/- 0.9 %	<+/- 0.8 %	< 0.6 % ¹	< 0.5 %	< 0.3 %

The performance values were determined using distilled water at a constant temperature (± 0.5 °C) between 20 and 25 °C in accordance with ISO 8655.

1) measured at 10 % nominal volume



Sarpette® M ordering information

Description	Volume range [µl]	Order no.	Packaging	Compatible SARSTEDT pipette tips
Sarpette® M single-channel	0.1 - 2 µl	90.3100.002	1 pc./case	70.3010.xxx 70.3020.xxx 70.3021.xxx 70.3010.xxx
	0.5 -10 µl	90.3100.010	1 pc./case	70.3020.xxx 70.3021.xxx 70.3020.xxx
	2 - 20 µl	90.3100.020	1 pc./case	70.3021.xxx 70.3030.xxx
	10 - 100 µl	90.3100.100	1 pc./case	70.3031.xxx 70.3030.xxx
	20 - 200 µl	90.3100.200	1 pc./case	70.3031.xxx 70.3050.xxx
	100 - 1,000 µl	90.3100.000	1 pc./case	70.3060.xxx 70.1183.102
	0.5 -5 ml	90.3100.555	1 pc./case	70.1183.002 70.1187.102
Sarpette® M 8-channel	1 - 10 ml	90.3100.111	1 pc./case	70.1187.002 70.3010.xxx
	0.5 -10 µl	90.3108.010	1 pc./case	70.3020.xxx 70.3021.xxx 70.3030.xxx
	20 - 200 µl	90.3108.200	1 pc./case	70.3031.xxx 70.3030.xxx
	40 - 300 µl	90.3108.300	1 pc./case	70.3031.xxx 70.3040.xxx
Sarpette® M 12-channel	0.5 -10 µl	90.3112.010	1 pc./case	70.3010.xxx 70.3020.xxx 70.3021.xxx
	20 - 200 µl	90.3112.200	1 pc./case	70.3030.xxx 70.3031.xxx
	40 - 300 µl	90.3112.300	1 pc./case	70.3030.xxx 70.3031.xxx 70.3040.xxx

Pipette holders for Sarpette® M ordering information

Description	Order no.	Packaging
Pipette holder, rotatable, for 6 SARSTEDT Sarpette® M single-channel pipettes up to 1,000 µl	95.1760.200	1 pc. / inner box 4 pcs. / case
Pipette holder, for 3 SARSTEDT Sarpette® M single- and multi-channel pipettes	95.1760.201	1 pc. / inner box 4 pcs. / case

Visit <https://refillrevolution.sarstedt.com/refill-revolution-usa/> for more information on the SARSTEDT pipette tip range.

SARSTEDT, Inc.

1025 St. James Church Road

P.O. Box 468

Newton, NC 28658 - 0468

Toll free: +1 800 257 5101

Fax: +1 828 465 4003

customerservice.us@sarstedt.com

www.sarstedt.com



SARSTEDT