

# S-Monovette® Serum Gel LightPROTECT

Light protection – from blood collection to archiving




- Light protection from the outset
- Light protection across the UV range
- More sustainability at the workplace by avoiding aluminium foil
- Usual convenient specimen identification as labels are always legible
- Process-compliant pre-analytics



**Pre-analytics**


The S-Monovette® Serum Gel LightPROTECT eliminates the need to wrap the blood sample in aluminium foil. Direct legibility and thus usual convenient specimen identification are maintained. This also facilitates an effortless contribution to more sustainability at the workplace by conserving resources (aluminium foil).

Light protection is guaranteed in the relevant UV range (190 nm – 380 nm). Average light permeability in this range is less than 1%.

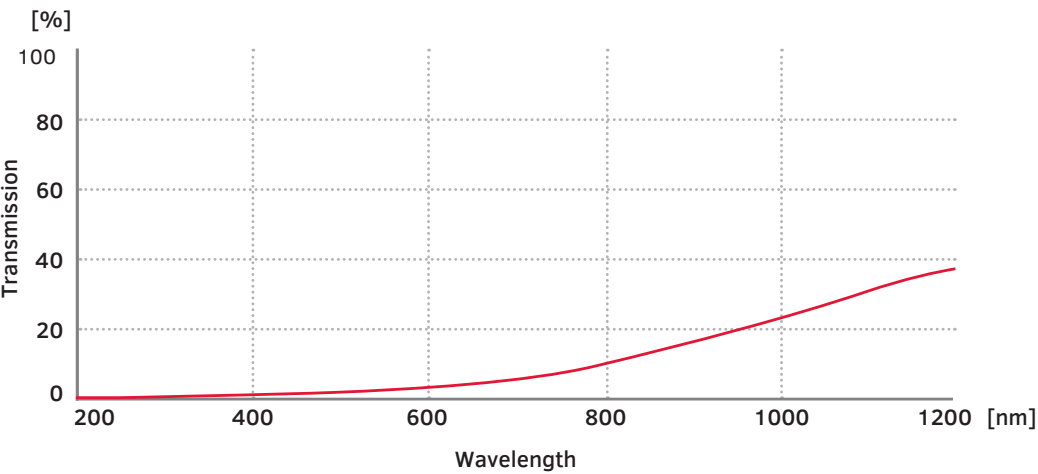

Transmission of the S-Monovette® Serum Gel LightPROTECT tube wall

### S-Monovette® Serum Gel LightPROTECT

Length / Ø	Colour	Case	Volume	Order no.
75 x 13 mm		500 / box 50 / inner box	4.0 ml	04.1925.003

### SARSTEDT AG & Co. KG

Sarstedtstraße 1  
D-51588 Nümbrecht  
Phone: +49 2293 305 0  
[export@sarstedt.com](mailto:export@sarstedt.com)  
[www.sarstedt.com](http://www.sarstedt.com)

