

*To make your
workflow excellent.*



Urine Drainage

Product solutions for the individual patient care



A06 W



DEHP free



- Needle-free and flat sampling site
- Attachment of the exchange bag with Bayonet-lock
- Universal attachment with elastic straps and additional hooks
- Exchange bags are available with 2,000 ml and 3,000 ml capacity




Precise Reading

The precise and easily readable graduations facilitate accurate output readings. The total measuring volume of 500 ml and the clear and concise graduations enable an accurate balancing of large urine volumes. Volumes exceeding 500 ml are automatically diverted through the integral overflow chamber and drained into the collection bag.

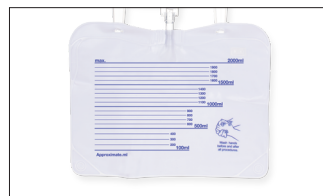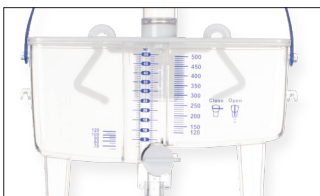

Hygienic & user-friendly exchange bag

The Bayonet-Lock enables a quick and hygienic exchange of the bag and a secure connection to the measuring chamber.

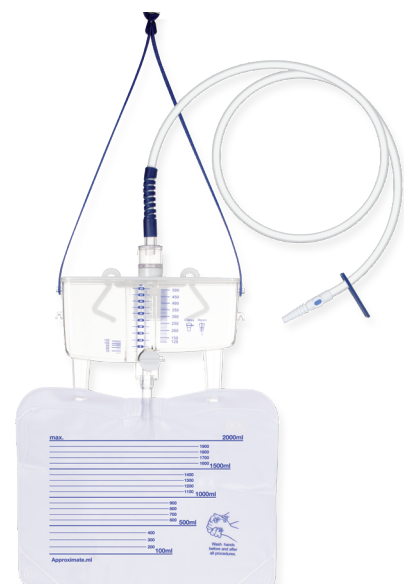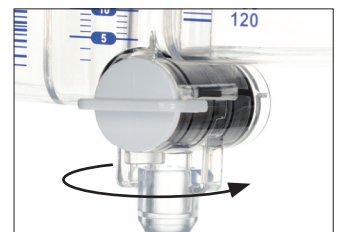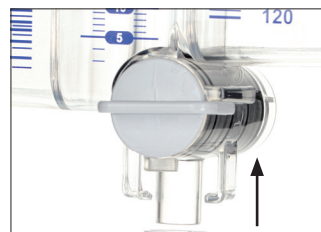

Needle-free sampling

Needle-free, flat sampling site for secure disinfection and sampling.



UD 3plus

DEHP free




- Needle-free and flat sampling site
- Large double hook
- 2 tube lengths: 100 cm and 120 cm
- 2,000 ml capacity

The System for the immobile patient

Closed drainage system for the hygienic care of catheterised patients with long-term confinement.

Large double hook

The large double hook allows easy and secure attachment to the bed rails.

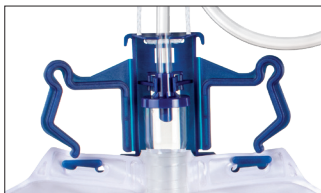

Drip Chamber

With hydrophobic ventilation for optimal pressure compensation and urine flow. The drip chamber enables optimal protection against ascending migration of bacteria. The horizontal membrane reflux valve ensures optimal urine flow from small to large volumes and prevents urine reflux regardless of the drainage bag position.



Needle-free sampling

Needle-free, flat sampling site for secure disinfection and sampling.



TUR-BAG

DEHP free



- Needle-free and flat sampling site
- Available with and without pumpball
- 4,000 ml capacity

The bag for bladder irrigation

A closed urine drainage system for the post-operative bladder irrigation.

Pumpball

With pump ball for controlled postoperative bladder irrigation, especially after prostate resection.



Needle-free sampling

Needle-free, flat sampling site for secure disinfection and sampling.



Drainage bag

DEHP free



- Universal step connector with integrated pierceable membrane for urine sampling
- Large diameter drainage tubing for optimal urine flow
- 2,000 ml capacity

Optimal Patientcare with the drainage bag

The drainage bag features an extra strong foil and a bacteria-proof, hydrophobic bag ventilation.

T-Tap drainage valve

Easy-to-handle T-Tap drainage valve and effective protection from contamination.



Needle-free sampling

Needle-free, flat sampling site for secure disinfection and sampling.



UroStretch



- No contact of the leg bag with the skin
- Full length fixation of leg bag surface
- Reduction of liquid noise

UroStretch leg bag attachment

The stretch fabric, tubular bandage for the attachment of urine leg-bags. Usable with all leg bags up to 750 ml. UroStretch is reusable, washable in warm water (40 °C), and offers high user comfort.



Ordering information

Urine Measuring System A06 W

Description	Pcs. / Box	Ref
Urine measuring system A06 W 500 ml with 2.0 l exchange bag without drainage tap, with elastic straps and additional hooks, individually packed, sterile	10	74.5681.200
Exchange bag 2.0 l for A06 W without drainage tap, individually packed, sterile	50	74.5684.001
Exchange bag 2.0 l for A06 W with drainage tap, individually packed, sterile	50	74.5684.002
Exchange bag 3.0 l for A06 W without drainage tap, individually packed, sterile	50	74.5685.001
Exchange bag 3.0 l for A06 W with drainage tap, individually packed, sterile	50	74.5685.002

UD 3plus

Description	Pcs. / Box	Ref
UD 3plus urine drainage system 2.0 l, with double hook, tube length 100 cm, individually packed, sterile	20	74.5330.500
UD 3plus urine drainage system 2.0 l, with double hook, tube length 120 cm, individually packed, sterile	20	74.5330.520

TUR-BAG

Description	Pcs. / Box	Ref
TUR-BAG 4.0 l with pump ball, tube length 100 cm, individually packed, sterile	10	74.5340.500
TUR-BAG 4.0 l without pump ball, tube length 100 cm, individually packed, sterile	20	74.5340.501

Drainagebeutel

Description	Pcs. / Box	Ref
Drainage bag 2.0 l, with anti-reflux valve, T-Tap drainage valve, tube length 100 cm, individually packed, sterile	50	74.5485.005

UroStretch

Description	Pcs. / Box	Ref
UroStretch, size: small	10	74.5767.001
UroStretch, size: medium	10	74.5767.002
UroStretch, size: large	10	74.5767.003
UroStretch, size: extra large	10	74.5767.004

If you have any questions,
we'll be happy to help!

Visit our website:
www.sarstedt.com

SARSTEDT AG & Co. KG

Sarstedtstraße 1
D-51588 Nümbrecht

Phone: +49 2293 305 0

export@sarstedt.com
www.sarstedt.com



Further information, videos and
sample orders:
sarstedt.com